Instructions for use

ConnectClip







Introduction to this booklet

This booklet guides you on how to use and maintain your new ConnectClip. Please read the booklet carefully, including the Warning section. This will help you get the most out of your new ConnectClip. If you have additional questions about the use of your ConnectClip, please contact your hearing care professional.

For your convenience, this booklet contains a color-coded navigation bar to help you navigate easily through the different sections.

Intended use

The ConnectClip is intended as a wireless receiver and transmitter of sound between hearing aids and the surrounding environment, including speech and wireless electronic audio devices. The ConnectClip also assists users of wireless hearing aids to control the functionality provided by their hearing solution.

The ConnectClip is intended to be used by children (> 36 months old) and adults together with their hearing solution.

IMPORTANT

Please familiarize yourself with the entire contents of this booklet before using your ConnectClip. It contains instructions and important information about the use and handling of your ConnectClip.

Introduction

The ConnectClip can be used both as a wireless headset for mobile phones and as a remote microphone/partner microphone.

20


1. Headset with mobile phone
Sounds from your phone
will be streamed to your
hearings aids through
the ConnectClip (A).
Your voice is picked

and sent to the mobile phone (B).

up by the ConnectClip

2. Remote/partner -microphone
The voice of a chosen speaker will

be transmitted wirelessly FROM the ConnectClip TO your hearing aids.

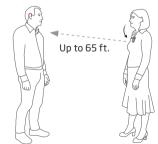

Contents

Pac	ckage contents	9
Ove	erview - main functions	10
Ove	erview - LED indicators	12
Get	ting started	14
Initial setup		
A. I	Pairing hearing aids with the ConnectClip	16
B. F	Pairing mobile phone with the ConnectClip	18
Daily use		
Dai	ly use: On/Off	20
Dai	ly use: Charging	21
Using Conne	ctClip with mobile phone	
Usi	ng ConnectClip with a mobile phone	22
Ens	suring ConnectClip is connected to the phone	23
Red	ceiving calls	24
Ma	king calls	26
Adv	vanced headset features	27


Listening to music	
Listening to music	28
■ Using ConnectClip as a remote microphone	
Using ConnectClip as a remote microphone	30
Selecting microphone mode	3.
Starting/stopping sound from remote microphone	37
Regulating level from remote microphone	33
Starting/stopping sound from remote microphone using your hearing aids	34
Cleaning the ConnectClip	35
■ Using ConnectClip as a remote control for hearing aids	
Using ConnectClip as a remote control for hearing aids	36
Other features and functions	37
Warnings	38
Technical information	44
International warranty	47


Package contents

The ConnectClip box consists of the following elements:



Overview - main functions

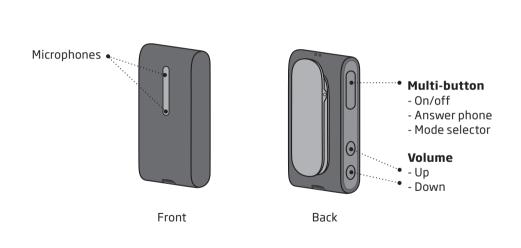

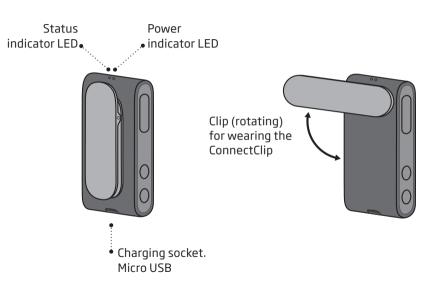

Overview - ConnectClip LED indicators


Status	Description	
•	Blinking yellow	Searching for hearing aids
•	Steady yellow	ldle. Hearing aids connected, Not connected to a phone
•	Blinking blue	Searching for a mobile phone (pairing)
•	Steady blue	Connected to a Mobile Phone
•	Steady green	In a call. Headset (HFP); Mobile
•	Blinking green	Incoming call
•	Steady orange	Music; Mobile Phone (A2DP)
•	Steady magenta	Remote microphone mode
•	Blinking magenta	Remote microphone reconnecting

Power	Description	
•	Steady red	Battery low
•	Steady green	Turned on - not charging
•	Blinking green	Turned on - charging
	Fading green	Turned on - charging complete
•	Blinking orange	Turned off - charging
•	Steady orange	Turned off - charging complete


Getting started

Before you can use the ConnectClip with your hearing aids, both the hearing aids and the ConnectClip need to be paired/connected. Using the ConnectClip with a mobile phone also requires the ConnectClip to be paired with/connected to your mobile phone.

1. Pairing hearing aid and ConnectClip

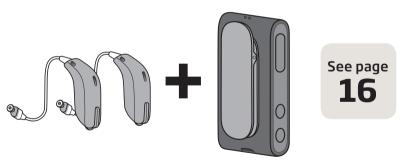

2. Pairing ConnectClip and mobile phone

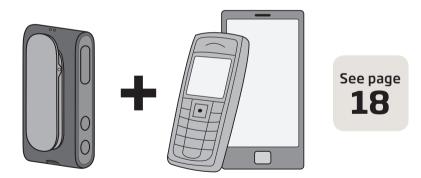


A. Pairing hearing aids with the ConnectClip


- Ensure ConnectClip is powered off (see page 20)
 Press and hold Multi-button for approx. 6 seconds until Power indicator OFF
- 2. Turn both hearing aids OFF (open battery drawer)
- 3. Turn both hearings aids ON (close battery drawer)



- **4. The hearing aids will now be ready** for first time connection with the ConnectClip. You now have 3 minutes to complete the pairing.
- 5. Turn ON ConnectClip Press and hold Multi-button for approx. 3 seconds until Power indicator shows STEADY GREEN

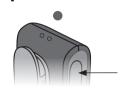




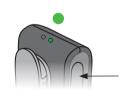

 ConnectClip will now search for your hearing aids.
 Status indicator BLINKING YELLOW




7. Place ConnectClip and both hearing aids close to each other (max distance 8 in.)

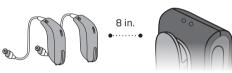

8. Connection established with hearing aid(s) is indicated by Status indicator STEADY YELLOW



B. Pairing mobile phone with the ConnectClip

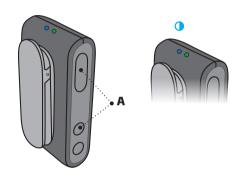

Place ConnectClip in pairing mode
Press and hold the Multi-button AND
Volume UP button for about 6 seconds
(A) until Status indicator BLINKS BLUE

The ConnectClip will be in pairing mode for about 3 minutes.



Enable pairing on mobile phone

Turn on Bluetooth. Search for and connect to/pair with ConnectClip.

This is normally in the settings/Bluetooth section of your phone.

Enter pin code


If your phone requests a pin code enter 0000 (four zeros).




Pairing complete

Pairing complete with mobile phone is indicated by status STEADY BLUE

Daily use: On/Off

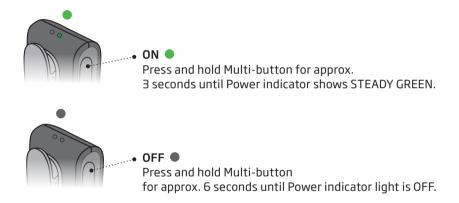


Daily use: Charging

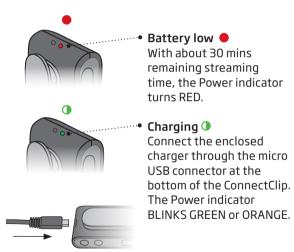

Charging indicators		
Charging (ConnectClip on)	Blinking	•
Charging complete (ConnectClip on)	Fading	
Charging (ConnectClip off)	Blinking	•
Charging complete (ConnectClip off)	Steady	•

20

Using ConnectClip with a mobile phone

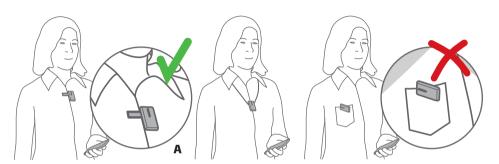

IMPORTANT

The ConnectClip should be worn as indicated for best possible performance. It is important that the microphone in the ConnectClip is pointing towards your mouth (A) in a vertical position, and that the ConnectClip is not covered by clothes etc. Maximum distance between ConnectClip and mobile phone is 32 ft.

Ensuring ConnectClip is connected to the phone



Select phone mode

Press and hold Multi-button for approx. 3 seconds to change mode. Mode change is confirmed by an audible indication in your hearing aids and ConnectClip Status indicator BLUE



Status	Mode
•	Idle Connected to hearing aids. Not connected to a mobile phone
•	Connected to mobile phone
•	Remote/partner microphone mode

Receiving calls



Phone connected

When the ConnectClip is turned on and is in range of your mobile phone it will automatically connect to the mobile phone.

The Status indicator is STEADY BLUE




Phone ringing

When the phone rings you will hear a ring indicator (and, if available, the ring-tone from your phone) in your hearing aids.

The Status indicator will BLINK GREEN ①




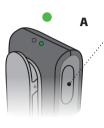

Answer call

Either click the Multi-button shortly (A) -or answer the call from your mobile phone.

Once the call is connected status will change to STEADY GREEN

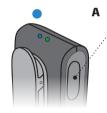

End call

Either click the Multi-button shortly (A) -or hang up the call from your mobile phone.

25

Status STEADY BLUE

Making calls







To make a call you only have to use the mobile phone. The call will automatically be connected to the ConnectClip.

Once the call is connected the status will change to STEADY GREEN •

Advanced headset features

Voice dialing / Voice command

When connected to a mobile phone, double press on the Multi-button activates voice dialing if supported by the phone.

Reject call

A long press on the Volume Down button, while the phone is ringing, rejects the call.

Mute hearing aid microphones

Press and hold Volume Down button for about 3 seconds to mute the hearing aid microphones.

Short press Volume Up or Down will unmute hearing aids

Listening to music



Wireless audio (A2DP)

If connected to a mobile phone or other device through Bluetooth, the music played on the phone/device will automatically be played in your hearing aids.


Status ORANGE

You do not need to operate anything on the ConnectClip to listen to music, but the ConnectClip must be turned on.

Adjusting volume


While listening to music, the level of the streamed signal can be adjusted using the Volume Up/Down buttons.

Short press on the Multi-button will stop sound. Use the controls on the phone to start the sound.

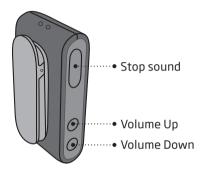

Status	Description
•	The ConnectClip is turned on - and streaming music from the phone

Using ConnectClip as a remote microphone

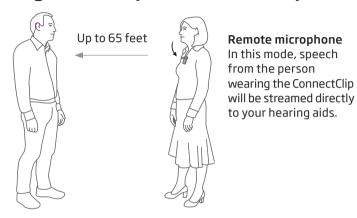

IMPORTANT

Ensure that the microphone is pointing to the speaker's mouth and that it is not covered by clothes etc.

Selecting microphone mode



• Select mic. mode

Press and hold Multi-button for approx. 3 seconds to change mode.


Mode change is confirmed by an audible indication in your hearing aids.

Status STEADY MAGENTA



Status	Description
•	The ConnectClip is turned on. In idle mode i.e. not connected to a phone
•	Connected to mobile phone
•	The ConnectClip is turned on - and in remote microphone mode

Starting/stopping sound from remote microphone

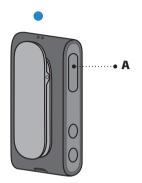

Stop sound from Microphone

Press and hold Multi-button (A) for approx. 3 seconds to stop the signal from the microphone to the hearing aids and switch the hearing aids back to their normal program.

The status inidcator will change to YELLOW or BLUE

32

Regulating level from remote microphone


Volume Up or Volume Down

The streamed signal can be adjusted using the Volume Up/Down buttons.

Mute Hearing aid microphones (hear only sound from the remote microphone).

Press and hold Volume Down for about 3 seconds to mute the hearing aid microphones.

Short press Volume Up or Down will unmute hearing aids.

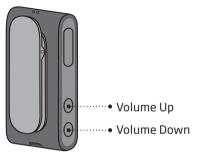

Starting/stopping sound from remote microphone using your hearing aids



P1 --> P2 --> TV --> Microphone

A long press on the hearing aid button will toggle to the next program in the hearing aids.

Once you reach the remote mirophone program you will hear the remote microphone program indicator.

To stop remote microphone sound simply toggle to back to the hearing aid program.

IMPORTANT

Please refer to your hearing aid's instructions for details on the functionality of your specific hearing aid solution.

Cleaning the ConnectClip

IMPORTANT

Use a soft, dry cloth to clean the ConnectClip. It must never be washed or immersed in water or other liquids.

Using ConnectClip as a remote control for hearing aids

When ConnectClip is NOT streaming sound, it doubles as a remote control for your hearing aids.



IMPORTANT

As remote control, maximum distance between ConnectClip & hearing aids is 6 ft.

36

Other features and functions

One connection at a time

When paired with several devices (up to eight), the ConnectClip will only connect to the first phone available/in range.

No calls when used as Remote microphone

When using the ConnectClip as Partner/Remote microphone - there will be no connection to the mobile phone.

Timeout (out of range) when in microphone mode

If the ConnectClip is out of range of the hearing aids - the hearing aids will return to their initial program after 90 seconds.

Clear Bluetooth pairings

Press and hold Multi-button, Volume Up AND Volume Down for 10 seconds to clear all Bluetooth pairings. Indicated by status BLINKING RED.

Marnings

You should familiarize yourself fully with the following general warnings and the entire contents of this booklet before using your ConnecClip to ensure personal safety and correct use. Consult your hearing care professional if you experience unexpected operations or events with your ConnectClip.

Choking hazards and risk of swallowing small parts

ConnectClip and its parts should be kept out of reach of children and anyone who might swallow these items or otherwise cause injury to themselves. If a small part is swallowed, see a doctor immediately.

Explosives

The power source in your ConnectClip has insufficient energy to cause fire in normal usage conditions. The ConnectClip has not been tested for compliance with international standards concerning explosive environments. We recommend you do not use your ConnectClip in areas where there is a danger of explosion.

Marnings

X-ray/CT/MR/PET scanning, electrotherapy and surgery

Remove your ConnectClip before X-ray/CT/MRI/PET scans, electrotherapy, surgery, etc. as your ConnectClip may be damaged when exposed to strong fields.

Heat and chemicals

The ConnectClip must never be exposed to extreme heat e.g. left inside a parked car in the sun. The ConnectClip must not be dried in microwave ovens or other ovens. The chemicals in cosmetics, hairspray, perfume, aftershave lotion, suntan lotion and insect repellent can damage the ConnectClip. Always remove your ConnectClip before applying such products and allow time to dry before use.

Interference

The ConnectClip has been thoroughly tested for interference according to the most stringent international standards. However, interference between the ConnectClip and other devices (e.g. some mobile telephones, citizens band devices, store alarm systems and other wireless devices) may occur. If this occurs, increase the distance between the ConnectClip and the interfering device.

Marnings

Use in traffic

Using the ConnectClip can distract you. Be aware of your surroundings while using the ConnectClip. Please refer to national and state traffic laws regarding the use of communication devices while in traffic.

Use on aircraft

Your ConnectClip contains Bluetooth. The ConnectClip should not be used on board an aircraft, unless specifically permitted by the flight personnel.

Active implants

Caution must be taken with active implants. If you have an active brain implant, please contact the manufacturer of your implanted device for information about the risk of disturbance.

Marnings

Rechargeable battery

Do not attempt to open the ConnectClip, as this may damage the battery and void the warranty.

Never attempt to replace the battery. If battery replacement is needed, please return your device to the supplier. The repair warranty is void if there are signs of tampering.

The safety of recharging batteries using a USB connector is determined by the external equipment. When the USB connector is connected to equipment plugged into a wall outlet, this equipment should be UL approved or comply with IEC-60065, IEC-60950 or other equivalent safety standards.

Connection to external equipment

The safety of the use ConnectClip when connected to external equipment (via auxiliary input cable, via USB cable, or directly), is determined by the external equipment. When connected to external equipment plugged into a wall outlet, this equipment must comply with IEC-60065, IEC-60950 or equivalent safety standards.

Marnings

Cross control

Your ConnectClip should be kept out of reach of children or anyone who might take control of the functionality, to ensure your personal safety and avoid incorrect use.

Use of lanyard

If by accident the lanyard is caught on something while wearing it, the lanyard is designed to break automatically. Therefore, never try to shorten or modify the lanyard in an unauthorized way. The lanyard should not be used by children under 36 months of age.

Possible side effects

The otherwise non-allergenic materials used in your ConnectClip may in rare cases cause a skin irritation or other side effects. Please consult a doctor if these conditions occur.

Marnings

Mobile phone compatibility

Please refer to www.oticon.com for a list of mobile phones that have been tested comprehensively by Oticon, and have been found to fulfil the minimum requirements for compatibility at the time of testing. However, Oticon takes no responsibility for such compatibility in general. Specifically, Oticon shall not be responsible for incompatibility resulting from future updates of the mobile phones, nor from events of incompatibility not covered by and/or found in these tests. Oticon takes no financial responsibility in terms of reimbursements of mobile phones that do not work as expected.

Before purchasing a new mobile phone, you are advised to consult the list mentioned above.

Technical information

Conditions of Use:

Operating humidity: 15% to 93% relative humidity, non-condensing.

Operating temperature: 5-40 °C / 41-104°F

Storage and transportation: -20-70 °C / -4-158°F

The ConnectClip contains a radio transceiver using Bluetooth EDR, Bluetooth Low Energy (BLE) and a proprietary short range radio technology all working at 2.4 GHz. The 2.4 GHz radio transmitter is weak and always below 17 dBm e.i.r.p. in total radiated power.

The ConnectClip complies with international standards concerning radio transmitters, electro-magnetic compatibility and human exposure.

Due to the limited space available on the ConnectClip, relevant approval markings can be found in this booklet.

USA and Canada

This device contains a radio module with the following certification ID numbers:

FCC ID: 2ACAHAC01 IC: 11936A-AC01

The device complies with Part 15 of the FCC Rules and with Industry Canada's licence-exempt RSSs.

Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This Class B digital apparatus complies with Canadian ICES-003.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not instal and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Redirect or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the manufacturer or an experienced radio/TV technician for help.

International warranty

Your ConnectClip is covered by an international limited warranty issued by the manufacturer for a period of 12 months from the date of delivery. This limited warranty covers manufacturing and material defects in the device itself, but not in accessories such as e.g. batteries, cables, power supplies, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, damage due to foreign objects entering the device or incorrect adjustments are NOT covered by the limited warranty and may void it. The above warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods. Your hearing care professional may have issued a warranty that goes beyond the clauses of this limited warranty. Please consult him/her for further information.

If you need service

Take your ConnectClip to your hearing care professional, who may be able to sort out minor problems and adjustments immediately.

The manufacturer declares that this ConnectClip (AC1A) is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU. The declaration of Conformity is available from the manufacturer.



Kongebakken 9 DK-2765 Smørum Denmark www.oticon.global













